### **Technical Sheet**

# GL PHOTOMETER 3.0 LS + FLICKER

High sensitivity and high sampling rate (125 kHz) photometer capable of extensive measurements of light flicker in addition to regular photometry. This device connects directly to PC by USB, to be controlled by GL Spectrosoft or dedicated API for external software use. Spectral response of the photometer precisely matches spectral sensitivity of human eye. This eliminates errors resulting from measurements of different type of light sources that are usual for this type of meters. LS version is prepared for goniometric measurements to work together with stray light eliminating tube. Such setup does not use optical elements responsible for cosine correction, what results in increased sensitivity.

#### Features:

- Spectral response of class A
- Increased sensitivity for goniophotometry
- 125 kHz 18 bit sampling rate
- Super wide dynamic range in a single setting mode
- Powered and controlled via USB connection




Application	High quality flicker measurements including long term sampling, goniometric measurements
MEASUREMENTS	
Measurement range	0.0001 lx 1 000 000 lx
Resolution	4 significant digits
Uncertainty of spectral response (f1')	Class A < 3 %
Sampling frequency	125 kHz
A/D conversion	18 bit
PROPERTIES	
Software	GL SPECTROSOFT or API for external software developers
Dimensions [H x W x D]	115 mm x 65 mm x 51 mm
Weight	350 g
Connectivity	USB A-B
Power	USB, 5V <200 mA
Tripod adapter	1/4"
ORDERING INFORMATION	
Case	$\checkmark$
USB cable	$\checkmark$
Part number:	No. 202294


Note: Instrument, firmware and software specification are subject to change without prior notice. All information included in GL OPTIC datasheets and product information available in any form are carefully prepared and included information believed to be true. Please note that discrepancies may occur due to text and/or other errors or changes in the available technology. We advise to contact GL Optic before the use of the product to obtain the latest product specification.





GL OPTIC POLSKA Sp. z o.o. ul. Poznańska 70, 62-040 Puszczykowo, Poland Phone: +48 61 819 40 03 | E-mail: office@gloptic.com www.gloptic.com **Technical Sheet** 


## GL PHOTOMETER 3.0 LS + FLICKER

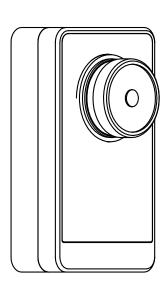

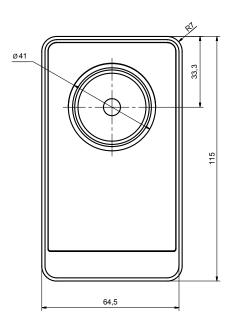

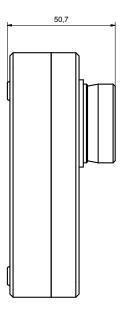

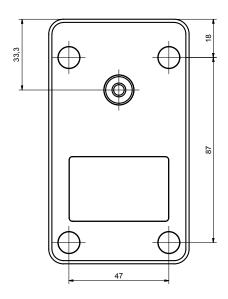

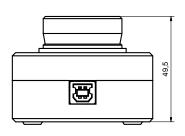



#### **GL OPTIC POLSKA Sp. z o.o.** ul. Poznańska 70, 62-040 Puszczykowo, Poland

ul. Poznańska 70, 62-040 Puszczykowo, Poland Phone: +48 61 819 40 03 | E-mail: office@gloptic.com www.gloptic.com